



# shouldfix light

Orthosis for shoulder  
immobilization

REF. 20413

Size  
Colour

universal size  
black

Fits right and left

## INDICATIONS

- Postoperative immobilisation following subacromial decompression, surgical removal of calcifications, sutures for the rotator cuff, revision surgery, arthrolysis, osteosynthesis, reverse shoulder replacement, and many other conditions
- All clinical cases that require early functional exercise (free movement of the elbow and wrist)
- Contusions
- Distortions

## CONTRAINDICATIONS

Currently no known



**Headquarter:** ORTHOSERVICE AG

Via Milano 7 · CH-6830 Chiasso (TI) · Switzerland  
Tel. 0041 (0) 91 822 00 88 · Fax 0041 (0) 91 822 00 89  
info@orthoservice.com · www.orthoservice.com

**Niederlassung Deutschland:** Orthoservice Deutschland GmbH

Flugstraße 8 · D-76532 Baden-Baden · Deutschland  
Tel. 0049 (0) 7221 991 39 11 · Fax 0049 (0) 7221 991 39 13  
info@orthoservice.de · www.orthoservice.de

**Sede italiana:** RO+TEN s.r.l.

**Sede legale:** Via Marco De Marchi, 7 · I-20121 Milano (MI) · Italia

**Sede operativa e amministrativa:**

Via Comasina, 111 · I-20843 Verano Brianza (MB) · Italia  
Tel. 0039 039 601 40 94 · Fax 0039 039 601 42 34  
info@roplusten.com · www.roplusten.com

Società soggetta a Direzione e Coordinamento (art. 2497bis CC):  
Orthoservice AG (CH) · 6830 Chiasso (TI) · Switzerland



UNI EN ISO 13485:2016  
UNI EN ISO 9001:2015



Organismo accreditato da ACCREDIA  
Body accredited by ACCREDIA




# shouldfix light

ORTHOSIS FOR SHOULDER  
IMMOBILIZATION

REF. 20413



ORTHOSERVICE ORTHOPEDIC PRODUCTS

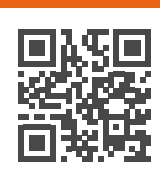

www.orthoservice.com







## Reduce pain and promote **healing**

**Shouldfix light** is used in the conservative and postoperative treatment of a range of shoulder conditions.

This brace has been designed to immobilise the shoulder in the most comfortable way possible yet without unduly restricting movement.

According to the latest theory about functional exercise, it has been shown that conventional immobilisation with internal rotation of the

arm can cause pain, tension and stiffness: this is why the basic configuration does not provide for fixing the forearm and wrist

**Shouldfix light** is the optimum solution in all clinical cases that require early functional exercise, such as distortion or contusion of the shoulder.



## ELBOW AND WRIST **MOBILITY**

The brace comprises an abdominal strap and an armband with Velcro®.

Made with nonslip inserts, the strap can be adapted to fit the patient and is easy to adjust. The armband is positioned above the elbow, and immobilises the arm while still allowing the elbow and wrist to move freely.

## An immobilised shoulder heals better with **freedom of movement**



## OPTIONAL **FOREARM RESTRICTION**

The patient's forearm and wrist can also be secured to the abdominal strap with easily positionable Velcro® bands.



The entire brace is made with a hassle-free fine, lightweight fabric that can even be worn against the skin without any problems, guaranteeing a comfortable fit and practical design. The brace can be easily washed.